

## Fast and Simple Free Fatty Acids Analysis Using UPC<sup>2</sup>/MS

Giorgis Isaac,<sup>1</sup> Michael D. Jones,<sup>1</sup> Besnik Bajrami,<sup>1</sup> Wassim Obeid,<sup>2</sup> James Langridge,<sup>3</sup> Patrick Hatcher<sup>2</sup>

<sup>1</sup>Waters Corporation, Milford, MA, USA

<sup>2</sup>Old Dominion University, Norfolk, VA, USA

<sup>3</sup>Waters Corporation, Manchester, UK

### APPLICATION BENEFITS

- Demonstrates the separation of free fatty acid (FFA) species based on chain length and number of double bonds
- No derivatization is required, which results in easier and fast sample preparation and eliminates artifact formation
- Organic phase lipid extract can be directly injected onto the system, saving time and reducing cost per analysis
- Less than three-minute chromatographic separation is up to 10X faster compared to GC/MS
- Unlike GC/MS, low volatile and very long chain fatty acids (>24 carbon atoms) can be easily analyzed with UPC<sup>2</sup>®

### WATERS SOLUTIONS

ACQUITY UPC<sup>2</sup>® System

TransOmics™ Informatics

Xevo® G2 QTof Mass Spectrometer

ACQUITY UPC<sup>2</sup> HSS Column

### KEY WORDS

Free fatty acids, UltraPerformance Convergence Chromatography™, UPC<sup>2</sup>, TransOmics, time-of-flight mass spectrometry, UPC<sup>2</sup>/MS/MS

### INTRODUCTION

Fatty acids, both free and as part of complex lipids, play a number of key roles in metabolism – as major metabolic fuel (storage and transport of energy), as essential components of all membranes, and as gene regulators. In addition, dietary lipids provide polyunsaturated fatty acids that are precursors of powerful locally acting metabolites, e.g., eicosanoids.

The common fatty acids of animal and plant origin have even-numbered chains of 16 to 24 carbon atoms with 0 to 6 double bonds. Nature provides countless exceptions, however, and odd- and even-numbered fatty acids with up to nearly 100 carbon atoms exist. In addition, double bonds can be of the *cis* (Z) and *trans* (E) configuration and there can be innumerable other structural features, including branch points, rings, oxygenated functions, and many more.

Fatty acid chains may contain one or more double bonds at specific positions (unsaturated and poly unsaturated with *cis* (Z) or *trans* (E) configuration) or they may be fully saturated. The LIPIDMAPS systematic nomenclature for fatty acids indicates the location of double bonds with reference to the carboxyl group with “Δ”.<sup>1</sup> Fatty acid structures also contain a methyl group at one end of the molecule (designated omega, ω) and a carboxyl group at the other end. The carbon atom next to the carboxyl group is called α carbon and the subsequent one the β carbon. The letter “n” is also often used instead of ω to indicate the position of the double bond closest to the methyl end.<sup>2</sup> Figure 1 outlines the structures of different straight chain fatty acids.

The isolation of free fatty acids (FFA) from biological materials is a complex task and precautions should be taken at all times to prevent or minimize the effects of hydrolyzing enzymes. After isolation, the typical chromatographic methods for analyzing fatty acids include gas chromatography/mass spectroscopy (GC/MS) and liquid chromatography-tandem mass spectrometry (LC/MS/MS). However, there are shortcomings associated with each of these methods.

For example, GC methods require derivatization of the fatty acids to hydrolyze and convert to methyl esters, which is time-consuming and risks re-arrangement of the fatty acids during derivatization, leaving doubt as to whether the esters formed are from FFA or intact complex lipids. Moreover, the GC/MS analysis of low volatile, very-long-chain fatty acids with high molecular weight (>C24) is a problem even after fatty acid methyl ester (FAME) derivatization.

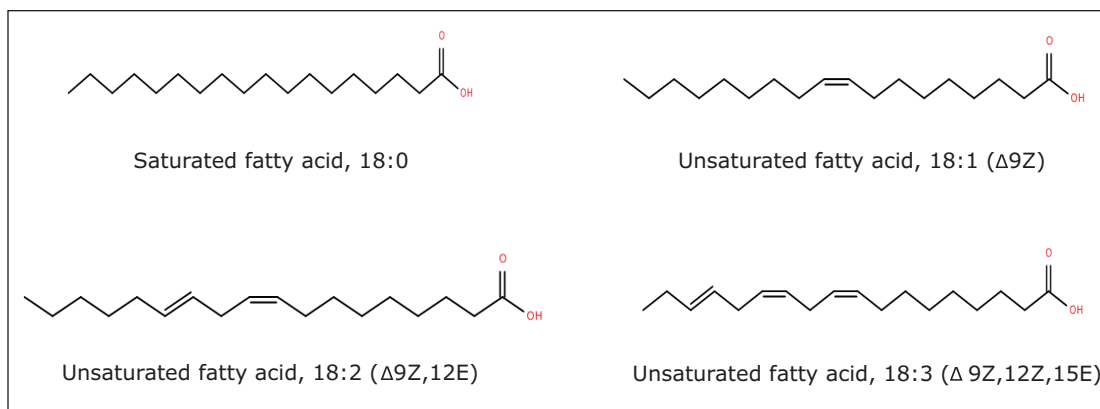


Figure 1. Structure and nomenclature of different straight chain fatty acids with a methyl and a carboxyl (acidic) end. Fatty acids may be named according to systematic or trivial nomenclature. One systematic way to describe the position of double bonds is in relation to the acidic end of the fatty acids; symbolized as  $\Delta$  (Greek delta) followed with numbers. All unsaturated fatty acids are shown with *cis* (Z) or *trans* (E) configuration of the double bonds.

In LC/MS methods, although no sample derivatization is required, the runs typically involve labor-intensive and time-consuming sample preparation, and utilize toxic organic solvents, which are expensive to purchase and dispose. In a typical reversed phase (RP) LC/MS analysis, the organic extracts containing all the lipids have to be evaporated and re-constituted in a more compatible injection solvent.

Thus, it would be beneficial to have streamlined methods for the separation and determination of fatty acids. Here, we present a rapid, high-throughput and efficient method for the separation and analysis of FFA using UltraPerformance Convergence Chromatography (UPCC, or UPC<sup>2</sup>) with mass spectrometry.

UPC<sup>2</sup> is a complementary, orthogonal separation technology that is taking its place alongside LC and GC. While all three use a stationary phase to interact with compounds of interest and a mobile phase to move compounds through the stationary phase and achieve separation, the techniques differ mainly by the mobile phases used.

GC is defined by using a gas as its mobile phase, LC is defined by using liquids as its mobile phase, and CC is defined by using both gas and liquids. It is this convergence of mobile phases in combination with a far greater choice of stationary phases that makes CC a powerful additional choice for laboratory scientists. Because UPC<sup>2</sup> can receive samples in organic solvents such as hexane and chloroform, it significantly simplifies the requirements for sample preparation, while maintaining all the advantages of RPLC.

Here, the analysis of fatty acids in the free form instead of FAME derivatives results in easier and faster sample preparation. The organic phase extract containing all the FFA can be injected directly into the system, which results in significant savings in sample preparation and analysis time, solvent costs, and solvent waste disposal. Additionally, artifact formation that can result from a derivatization procedure is eliminated.

**EXPERIMENTAL****Method Conditions****UPC<sup>2</sup> conditions**

System:	ACQUITY UPC <sup>2</sup>
Columns:	ACQUITY UPC <sup>2</sup> HSS C <sub>18</sub> SB 1.8 μm, 2.1 x 150 mm
Column temp.:	50 °C
Sample vial:	Total Recovery Vial (p/n 186000385C)
Sample temp.:	10 °C
Injection volume:	0.5 μL
Flow rate:	0.6 mL/min
Mobile phase A:	CO <sub>2</sub>
Mobile phase B:	Methanol in 0.1% formic acid
Make up:	Methanol in 0.1% NH <sub>4</sub> OH (0.2 mL/min)
Splitter:	Upchurch cross 1/16 PEEK

**Gradient**

Time (min)	%A (CO <sub>2</sub> )	%B	Curve
0.0	95	5	Initial
5.0	75	25	6
5.1	50	50	1
6.0	50	50	11
8.0	95	5	1

**MS conditions**

Mass spectrometer:	Xevo G2 QTof
Ionization mode:	ESI negative
Capillary voltage:	1.0 kV
Cone voltage:	30 V
Source temp.:	100 °C
Desolvation temp.:	500 °C
Cone gas flow:	10 L/h
Desolvation gas flow:	600 L/h
Acquisition range:	50 to 600 <i>m/z</i>

**Sample preparation****FFA standard mixtures**

Individual saturated FFA standards containing even carbon number C<sub>8</sub> to C<sub>24</sub> were purchased from Sigma. A complex model mixture of different FFA standards (GLC-85 in FFA form) was purchased from Nu-Chek Prep (Elysian, MN, USA). The list of FFA standards analyzed and other detailed information is provided in Table 1. A 1 mg/mL stock solution was prepared in chloroform, and 0.1 mg/mL working lipid mixtures were prepared in chloroform, then injected onto the UPC<sup>2</sup>/MS system.

**Algae and algaenan produced oils**

Oil produced from hydrous pyrolysis of algae and algaenan at low and high pyrolysis temperature were provided from Old Dominion University (Norfolk, VA, USA). Algae 1 and algaenan 1 were treated at a pyrolysis temperature of (310 °C); and Algae 2 and algaenan 2 were treated at a pyrolysis temperature of (360 °C).

Extraction of algaenan was performed by a modified extraction procedure. Briefly, lipids were removed from the algae by Soxhlet extraction with 1:1 (v/v) benzene/methanol solvent mixture for 24 hours. The residue was treated with 2N sodium hydroxide at 60 °C for two hours. The remaining residue was then washed excessively with deionized water, followed by treatment with Dowex 50W-x8 cation exchange resin to exchange any residual sodium. Finally, the solid was rinsed with deionized water. The oil samples were diluted 10 times in dichloromethane, and 1 μL was injected onto the UPC<sup>2</sup>/MS system.

**Data acquisition and processing**

When using multivariate data analysis for sample comparison, it is crucial that each sample is randomized and injected a minimum of three times to ensure that the data analysis is statistically valid. For this study, five replicates of each algae and algaenan oil extracts were acquired in MS<sup>E</sup> mode, an unbiased ToF acquisition method in which the mass spectrometer switches between low and elevated collision energy on alternate scans. Data analysis and FFA identification were performed using TransOmics Informatics for Metabolomics and Lipidomics (TOIML).

Compound	Formula	Neutral Mass	[M-H] <sup>-</sup>	Retention time (min)	Common Name	Description
1	C <sub>4</sub> H <sub>8</sub> O <sub>2</sub>	88.052429	87.045153	0.89	Butyric acid	C4:0
2	C <sub>6</sub> H <sub>12</sub> O <sub>2</sub>	116.083730	115.076453	0.96	Caproic acid	C6:0
3	C <sub>8</sub> H <sub>16</sub> O <sub>2</sub>	144.115030	143.107753	1.06	Caprylic acid	C8:0
4	C <sub>10</sub> H <sub>20</sub> O <sub>2</sub>	172.146330	171.139053	1.17	Capric acid	C10:0
5	C <sub>11</sub> H <sub>22</sub> O <sub>2</sub>	186.161980	185.154704	1.23	Undecylic acid	C11:0
6	C <sub>12</sub> H <sub>24</sub> O <sub>2</sub>	200.177630	199.170354	1.31	Lauric acid	C12:0
7	C <sub>13</sub> H <sub>26</sub> O <sub>2</sub>	214.193280	213.186004	1.41	Tridecylic acid	C13:0
8	C <sub>14</sub> H <sub>28</sub> O <sub>2</sub>	228.208930	227.201654	1.54	Myristic acid	C14:0
9	C <sub>15</sub> H <sub>30</sub> O <sub>2</sub>	242.224580	241.217304	1.67	Pentadecylic acid	C15:0
10	C <sub>16</sub> H <sub>32</sub> O <sub>2</sub>	256.240230	255.232954	1.80	Palmitic acid	C16:0
11	C <sub>17</sub> H <sub>34</sub> O <sub>2</sub>	270.255880	269.248604	1.97	Margaric acid	C17:0
12	C <sub>18</sub> H <sub>36</sub> O <sub>2</sub>	284.271530	283.264254	2.11	Stearic acid	C18:0
13	C <sub>20</sub> H <sub>40</sub> O <sub>2</sub>	312.302831	311.295554	2.41	Arachidic acid	C20:0
14	C <sub>22</sub> H <sub>44</sub> O <sub>2</sub>	340.334131	339.326854	2.70	Behenic acid	C22:0
15	C <sub>14</sub> H <sub>26</sub> O <sub>2</sub>	226.193280	225.186004	1.45	Physeteric acid	C14:1
16	C <sub>15</sub> H <sub>28</sub> O <sub>2</sub>	240.208930	239.201654	1.57		C15:1
17	C <sub>16</sub> H <sub>30</sub> O <sub>2</sub>	254.224580	253.217304	1.67	Palmitoleic acid	16:1
18	C <sub>17</sub> H <sub>32</sub> O <sub>2</sub>	268.240230	267.232954	1.81	10-HEPTADECENOIC Acid	C17:1 (Δ10)
19	C <sub>18</sub> H <sub>30</sub> O <sub>2</sub>	278.224580	277.217304	1.76	Gamma Linolenic Acid	C18:3 (Δ6,9,12)
20	C <sub>18</sub> H <sub>30</sub> O <sub>2</sub>	278.224580	277.217304	1.86	Linolenic Acid	C18:3 (Δ9,12,15)
21	C <sub>18</sub> H <sub>30</sub> O <sub>2</sub>	280.240230	279.232954	1.88	Linoleic Acid	C18:2
22	C <sub>18</sub> H <sub>34</sub> O <sub>2</sub>	282.255880	281.248604	1.98	Oleic Acid	C18:1
23	C <sub>18</sub> H <sub>34</sub> O <sub>2</sub>	282.255880	281.248604	1.98	Elaidic Acid	C18:1T
24	C <sub>20</sub> H <sub>32</sub> O <sub>2</sub>	304.240230	303.232954	1.93	Arachidonic acid	C20:4
25	C <sub>20</sub> H <sub>34</sub> O <sub>2</sub>	306.255880	305.248604	2.04	HOMOGAMMA LINOLENIC Acid	C20:3 (Δ8,11,14)
26	C <sub>20</sub> H <sub>34</sub> O <sub>2</sub>	306.255880	305.248604	2.14	11-14-17-EICOSATRIENOIC Acid	C20:3 (Δ11,14,17)
27	C <sub>20</sub> H <sub>36</sub> O <sub>2</sub>	308.271530	307.264254	2.17	11-14-EICOSADIENOIC Acid	C20:2 (Δ11, 14)
28	C <sub>20</sub> H <sub>38</sub> O <sub>2</sub>	310.287180	309.279904	2.24	11-EICOSENOIC Acid	C20:1 (Δ11)
29	C <sub>22</sub> H <sub>32</sub> O <sub>2</sub>	328.240230	327.232954	2.09	Docosahexaenoic Acid	C22:6
30	C <sub>22</sub> H <sub>40</sub> O <sub>2</sub>	336.302831	335.295554	2.46	Docosadienoic Acid	C22:2
31	C <sub>22</sub> H <sub>38</sub> O <sub>2</sub>	338.318481	337.311204	2.54	Erucic Acid	C22:1
32	C <sub>24</sub> H <sub>46</sub> O <sub>2</sub>	366.349781	365.342504	2.83	Nervonic acid	C24:1

Table 1. A list of analyzed saturated and unsaturated standard FFA mixtures with corresponding retention time determined from Figure 3A.

## RESULTS AND DISCUSSION

### Analysis of saturated FFA standards

Figure 2 shows the separation of saturated FFA with carbon chain length  $C_8$  to  $C_{24}$ . The ACQUITY UPC<sup>2</sup> High Strength Silica (HSS)  $C_{18}$  SB 1.8  $\mu$ m, 2.1 x 150 mm Column provides an RP-like separation that results in effective separation of the different FFA species. The gradient is run under acidic conditions using a small percentage of formic acid (0.1% v/v in methanol) to improve the peak shape and decrease peak tailing.

The ACQUITY UPC<sup>2</sup> method is 10X faster (only a three-minute run) than GC/MS and RPLC methods, and uses less toxic and cheaper  $CO_2$  as a solvent. A typical lipidomics study involves the analysis of thousands of biological samples, and the additional speed allows for large sample sets to be analyzed efficiently, improving the overall power of the experiment.

The FFA lipid molecular species separation mechanism is mainly based on hydrophobic interaction of the FFA carbon numbers and number of double bonds with the HSS  $C_{18}$  SB material. Therefore, the elution order of the FFA species depends on the length and the number of double bonds on the fatty acid chain. Thus, the longer and the more saturated the acyl chain length the longer the retention time.

The co-solvent mobile phase B (methanol in 0.1% formic acid) can be optimized to increase the chromatographic resolution and peak capacity. The higher the percentage of the co-solvent, the shorter the retention time and the narrower the peaks. However, when analyzing a complex biological sample containing saturated and unsaturated FFA species with different carbon chain length, peak capacity is important in order to reduce coeluting lipid species. The co-solvent gradient 5% to 25% methanol in 0.1% formic acid was used for further analysis.

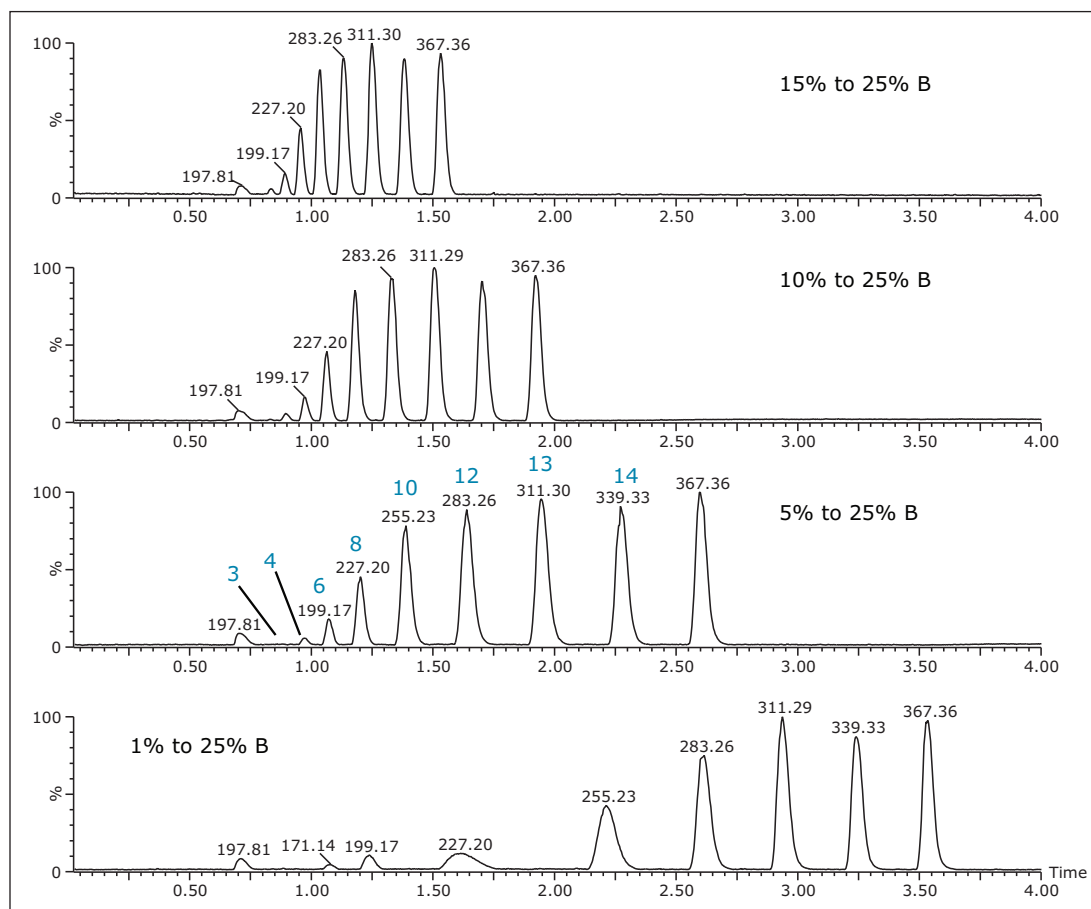


Figure 2. The separation of saturated FFA with carbon chain length  $C_8$ - $C_{24}$  with various co-solvent gradient. For the lipid ID, see Table 1.

## Analysis of complex saturated and unsaturated FFA standards GLC-85

Reversed-phase chromatography separates lipids according to both chain-length and degree of unsaturation. The problem lies in the fact that the dual nature of the reversed-phase separation process (a double bond in the fatty acyl chain reduces the retention time and the fatty acyl chain length increases the retention time) can hamper the analysis of real samples; the number of components is often so great that identification becomes difficult due to coelution (Figures 3A and B).

On the other hand, by using the precursor exact mass, corresponding product ion information and ion mobility (separation of lipid ions in the gas phase according to their size and molecular shape), each coeluting peak can be extracted and identified.

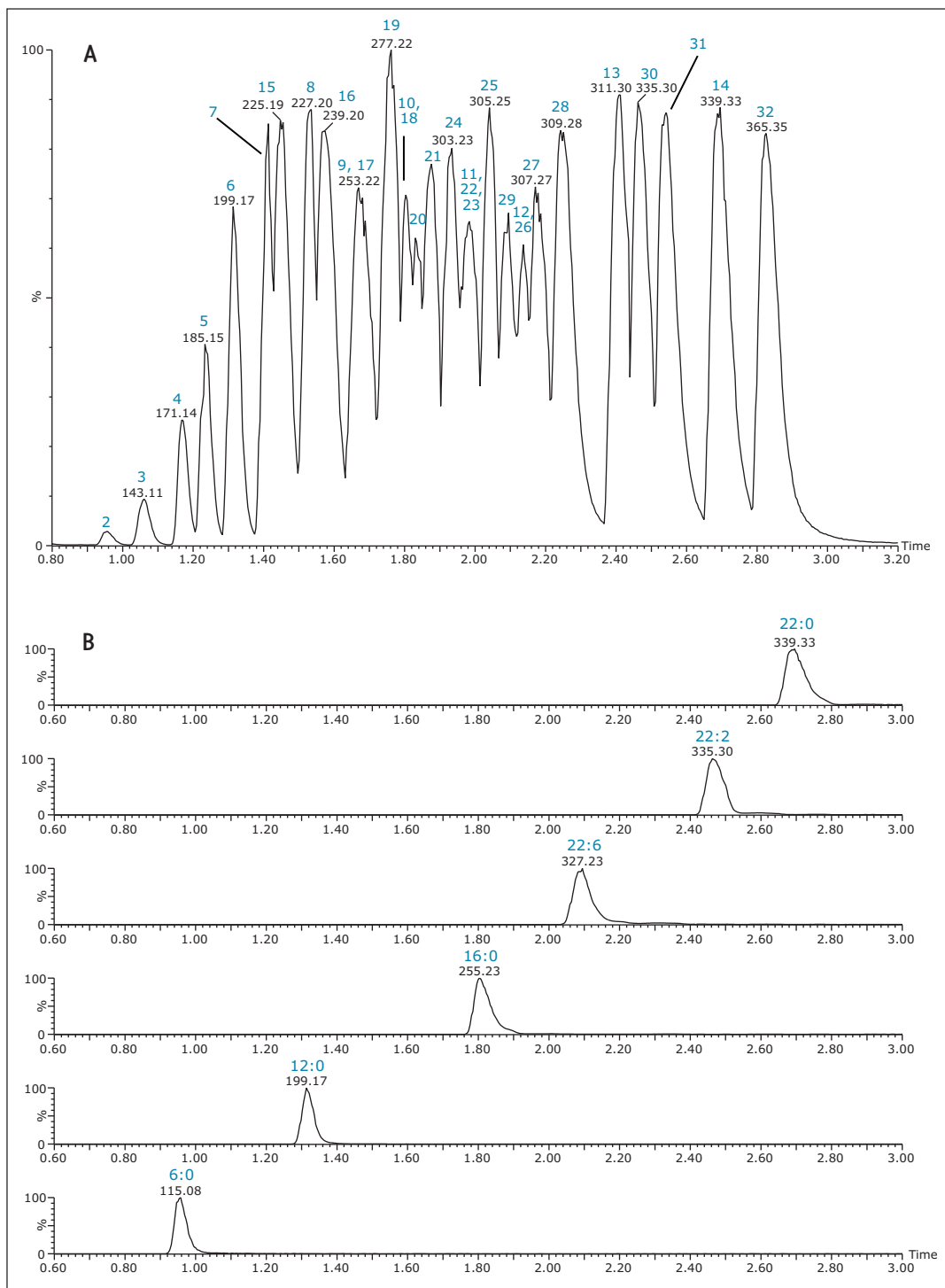


Figure 3. A) The separation of complex standard mixture that contains saturated, unsaturated, short and long chain 32 different FFA species. B) The separation depends on both chain length and degree of unsaturation. In an RP separation, the fatty acyl chain length increases the retention time and the number of double bonds in the fatty acyl chain decreases the retention time. For the lipid ID, see Table 1.

Another benefit of the method is the ability to separate between lipid isomers. FFA can have different biological functions based on the double bond position (e.g., omega-3 and omega-6). Figure 4 shows the separation of FFA isomers based on the position of the double bond. The separation of 18:3 ( $\Delta 6,9,12$ ) and 18:3 ( $\Delta 9,12,15$ ); and 20:3 ( $\Delta 8,11,14$ ) and 20:3 ( $\Delta 11,14,17$ ) isomers have been observed.

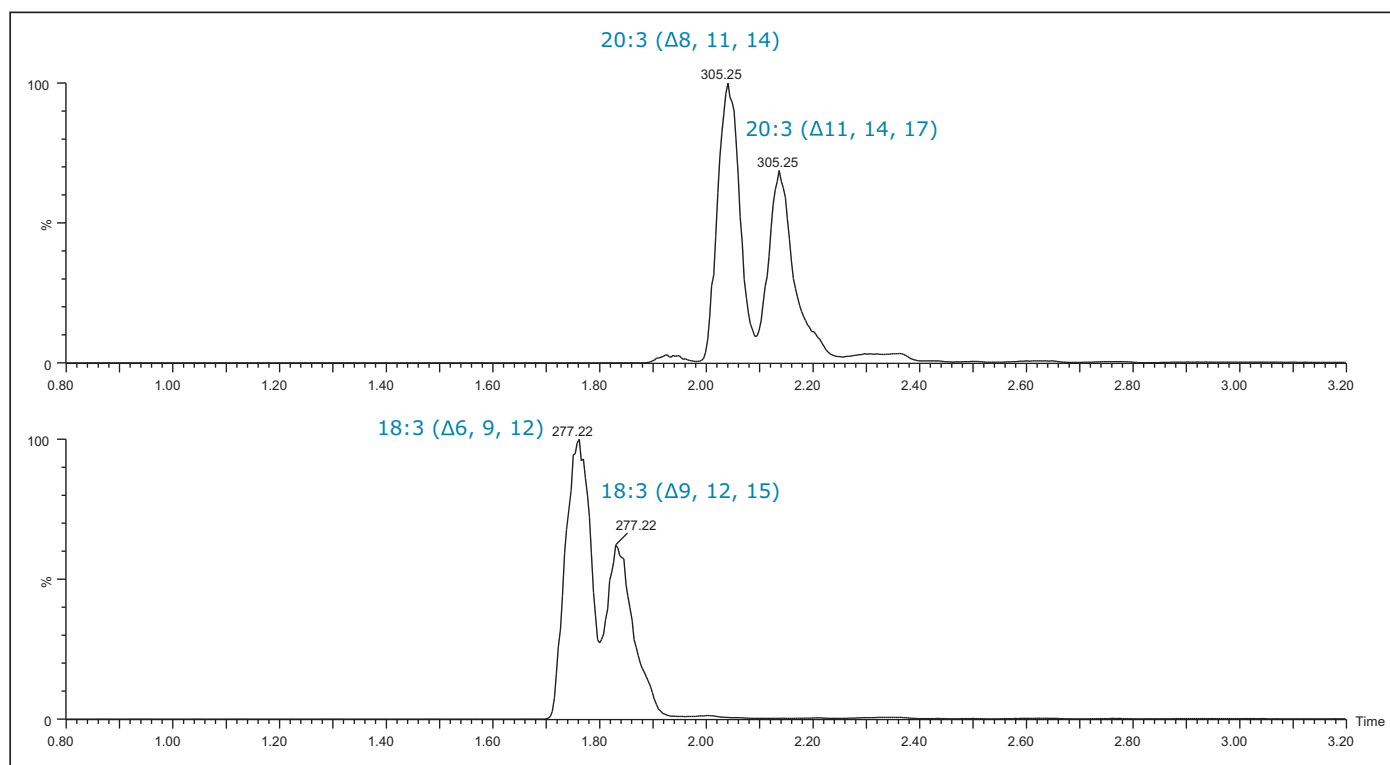


Figure 4. Extracted ion chromatogram (from figure 3) showing the separation of isobaric lipid species based on the position of the double bond.

## Biological application and data analysis using TransOmics

The developed UPC<sup>2</sup>/Xevo G2 QTof MS method was applied with minor modifications for the profile of FFA in algae and algaenan extracts treated at low (310 °C) and high (360 °C) pyrolysis temperatures.

Algaenan is a non-hydrolyzable, insoluble biopolymer in the cell walls of several green freshwater and marine microalgae.<sup>3</sup> Figure 5 shows a representative chromatogram from algaenan 1 with the UPC<sup>2</sup> conditions used for the analysis. For complete analysis of the data, set the gradient 1% to 10% co-solvent mobile phase B (methanol in 0.1% formic acid) in 10 minutes was used.

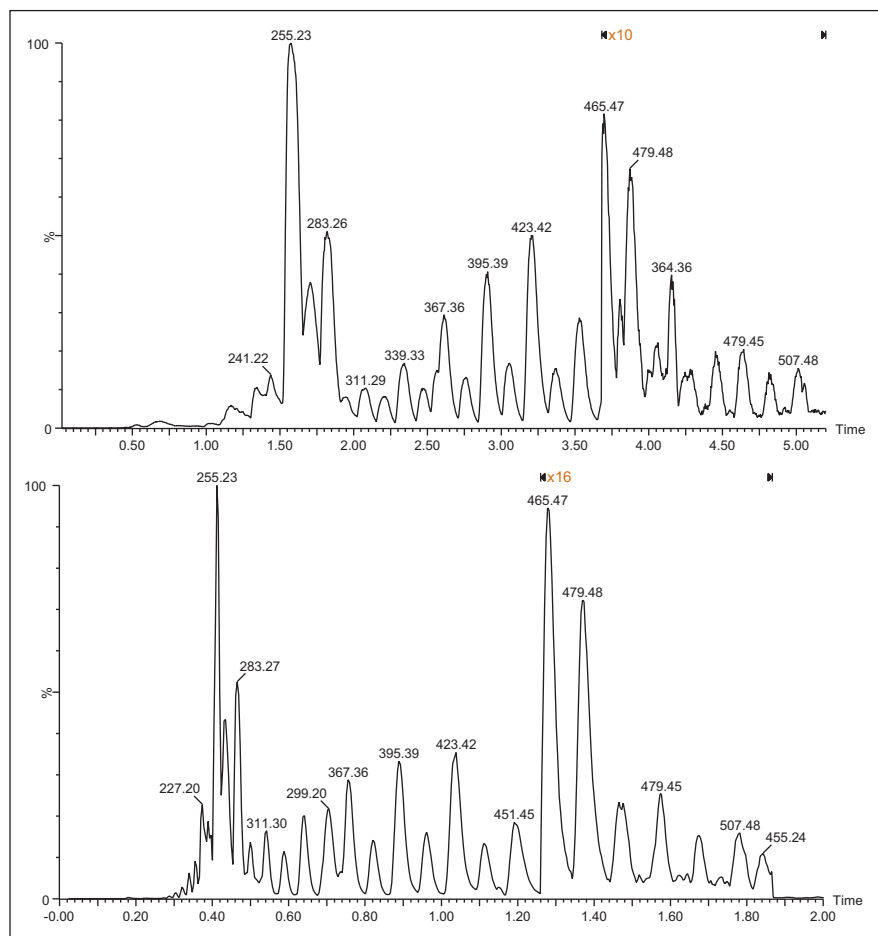


Figure 5. Representative chromatogram from algaenan 1 with various co-solvent gradients (top 1% to 10% methanol in 10 minutes, lower 5% to 20% methanol in 10 minutes). (UPC<sup>2</sup> conditions: HSS C<sub>18</sub> SB (2.1 x 100 mm), flow rate= 1.5 mL/min. The other UPC<sup>2</sup> conditions are described in the method conditions).



The lipid profiles of the algae and algaenan oil were investigated using TransOmics (TOIML) Software to determine the pattern and composition of FFA at two different pyrolysis temperatures. Differential analysis of results across different treatments can quickly be performed, thereby facilitating identification and quantitation of potential biomarkers. The software adopts an intuitive workflow approach to performing comparative UPC<sup>2</sup>/Xevo G2 QTof MS metabolomics and lipidomics data analysis.

The workflow starts with UPC<sup>2</sup>/MS raw data file loading, then retention time alignment and deconvolution, followed by analysis that creates a list of features. The features are then identified with compound searches and explored using multivariate statistical methods.

Principal component analysis (PCA) was used in the first instance to identify the combination of the FFA species that best describe the maximum variance between algae 1, algae 2, algaenan 1, and algaenan 2 oils (Figure 6). The PCA plot showed excellent technical UPC<sup>2</sup>/MS measurements. The PCA plot effectively displays the inter-sample relationships in multi-dimensional hyperspace, with more similar samples clustering together and dissimilar samples separated.<sup>4</sup>

The clustering in Figure 6 indicates that algae 1 and algaenan 1 are different, but algae 2 and algaenan 2 have more similarity in their FFA compositions after high pyrolysis temperature treatment. Orthogonal projections latent structure discriminant analysis (OPLS-DA) binary comparison can be performed between the different sample groups (algae 1 vs. algae 2, algaenan 1 vs. algaenan 2, algae 1 vs. algaenan 1, and algae 2 vs. algaenan 2) to find out the features that change between the two groups.

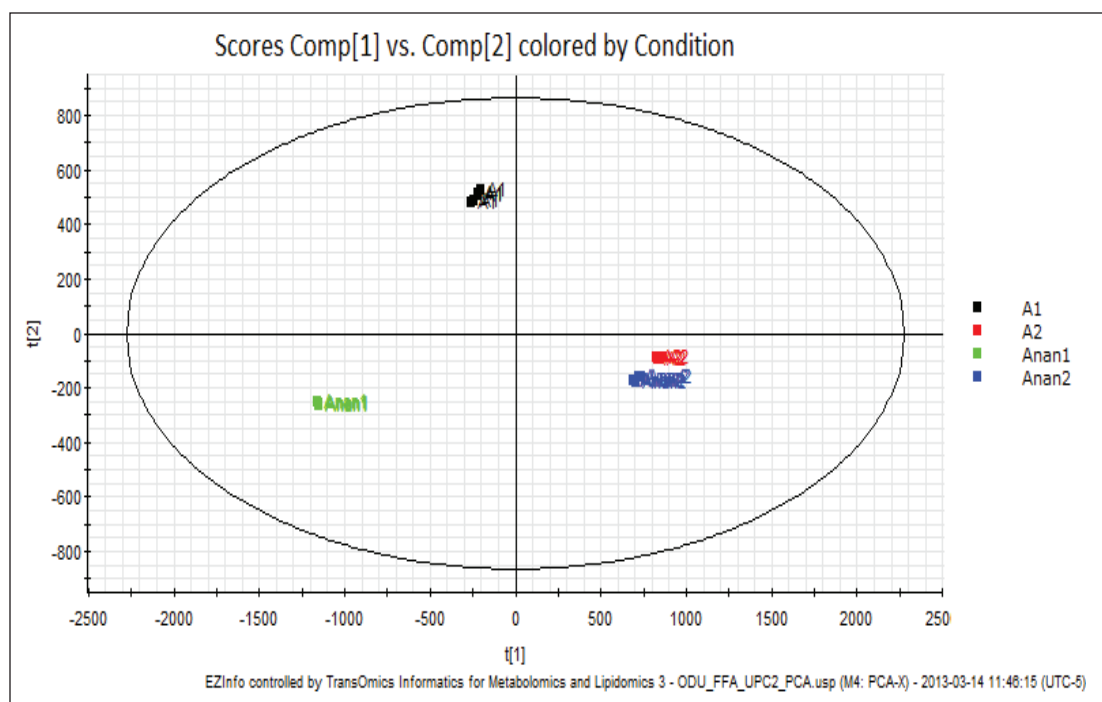


Figure 6. Principal component analysis of algae and algaenan oil extracts treated at low and high pyrolysis temperature. (A1= algae at low pyrolysis temperature A2= algae at high pyrolysis temperature Anan1= algaenan at low pyrolysis temperature Anan2= algaenan at high pyrolysis temperature).

As an example, the OPLS-DA binary comparison between algae 1 vs. algae 2 is shown in Figure 7A. As shown in the S-plot, the features that contribute most to the variance between the two groups are those farthest from the origin of the plot, highlighted in red (Figure 7B). These selected features can be exported to TransOmics for further identification. This helps the researcher focus on the features/compounds that change between samples instead of spending time on the whole data set.

Figures 7C and 7D show representative trend plots that change most between algae 1 and algae 2. Figure 8A shows the ion map, mass spectrum, and chromatogram across all the runs for FFA 29:0. This view allows to review compound measurements such as peak picking and alignment to ensure they are valid across all the runs. Figure 8B shows the normalized abundance of FFA 29:0 across all the conditions. FFA 29:0 is elevated in algae 1 compared to algae 1, algae 2, and algae 2; however, there is no significant difference between algae 2 and algae 2. Detailed investigation and comparison between algae 1 and algae 2 showed that algae 1 contains elevated levels of short (C9:0 to C13:0) and long (C31:0 to C37:0) chain FFA, whereas algae 2 contains elevated levels of medium (C14:0-C29:0) chain FFA. Similarly, the comparison between algae 1 and algae 2 showed that algae 1 contains elevated levels of long (C28:0 to C37:0) chain FFA, whereas algae 2 contains elevated levels of short and medium (C9:0 to C27:0) chain FFA.

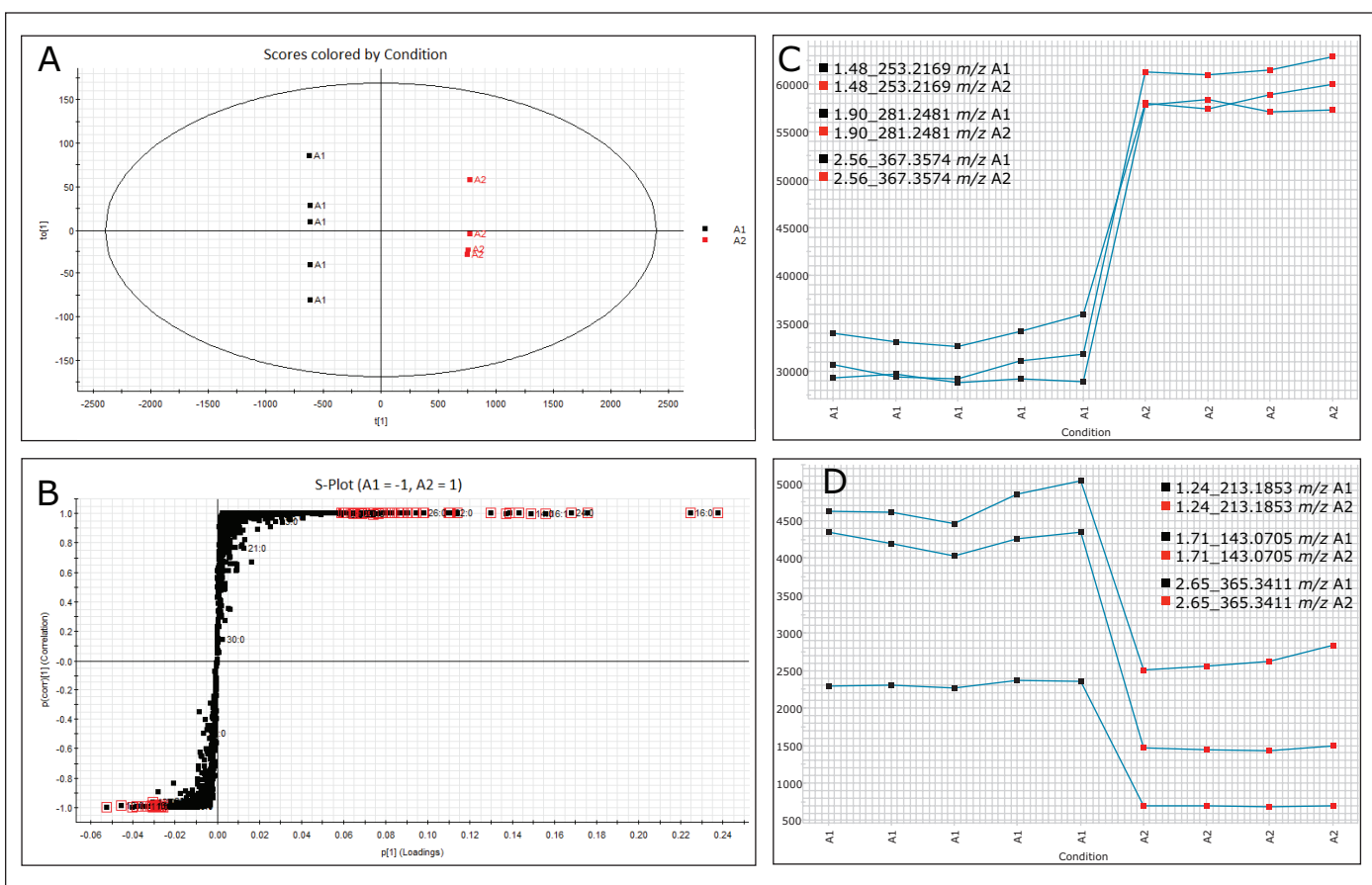


Figure 7. (A) OPLS DA plot between algae 1 and algae 2 group difference. (B) S-plot indicating the major features (highlighted in red) that contribute to the group difference between algae 1 and algae 2. (C) Representative trend plot showing the major up-regulated 16:1, 18:1, and 24:0 FFA in A1 (D) Representative trend plot showing the major up-regulated 8:0, 13:0, and 24:1 FFA in A2. (A1= algae at low pyrolysis temperature A2= algae at high pyrolysis temperature).

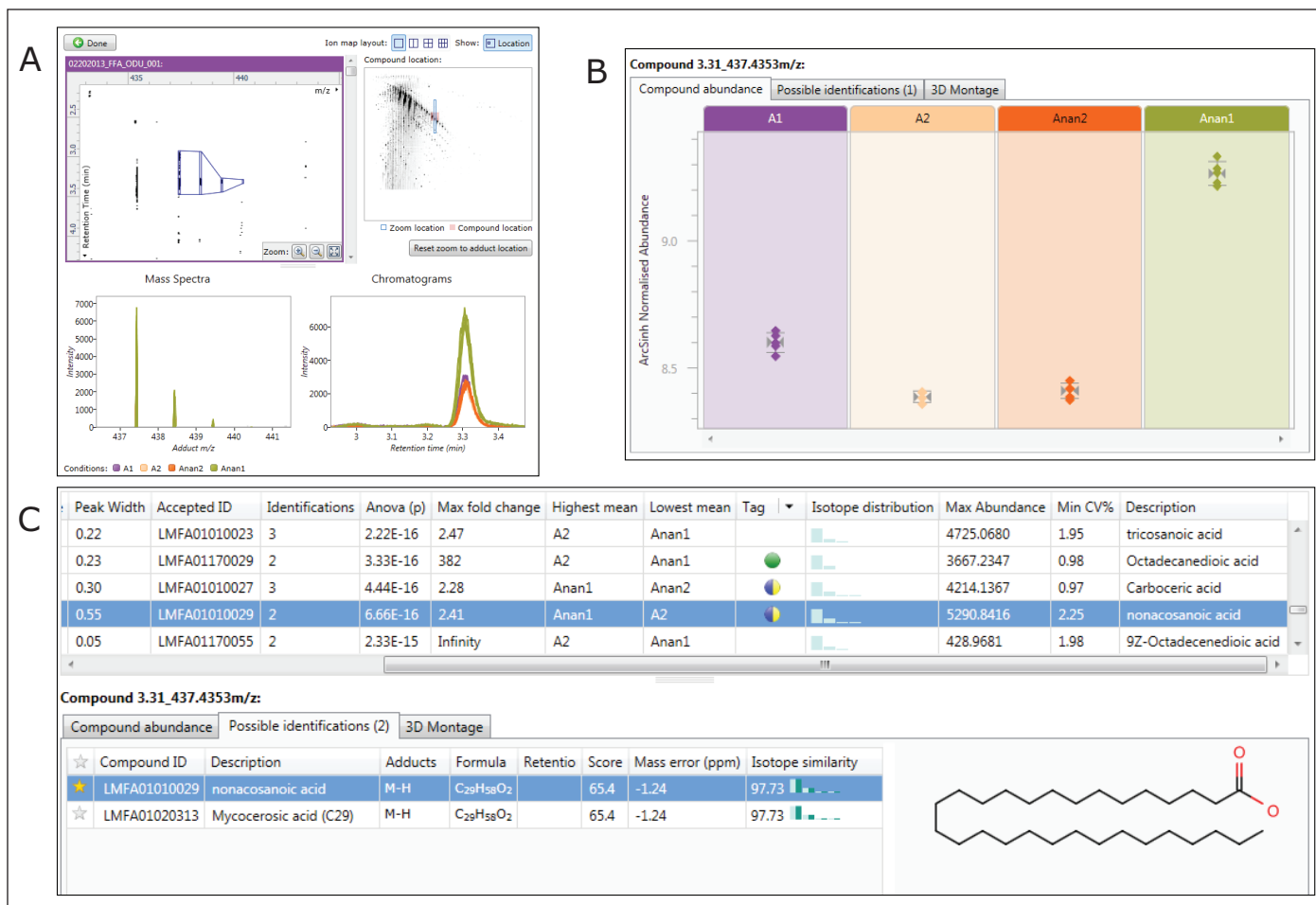


Figure 8. (A) Selected FFA 29:0 showing its ion map, mass spectrum, and chromatogram across all the runs. (B) Normalized abundance of FFA 29:0 across all the conditions. (C) Identification can be performed by means of local or web-based database search. In this example, the feature with retention time and exact mass pair 3.31\_437.4353 is identified as nonacosanoic acid (29:0 FFA). (A1= algae at low pyrolysis temperature, A2= algae at high pyrolysis temperature; Anan1= algaenan at low pyrolysis temperature Anan2= algaenan at high pyrolysis temperature).

Identification can be performed by means of local or web-based (such as LIPID MAPS, HMDB, and METLIN) compound searches based on retention time, low energy exact mass, high energy fragment ion, theoretical isotope pattern distribution, and collision cross section area (CCS) (Figure 8C). In this example, the feature with retention time and exact mass 3.31\_437.4353 is identified as nonacosanoic acid (29:0 FFA) based on retention time, low energy exact mass, and theoretical isotope pattern distribution. Figure 9 shows the expression and abundance profile of selected features according to their relative similarity between the different groups.

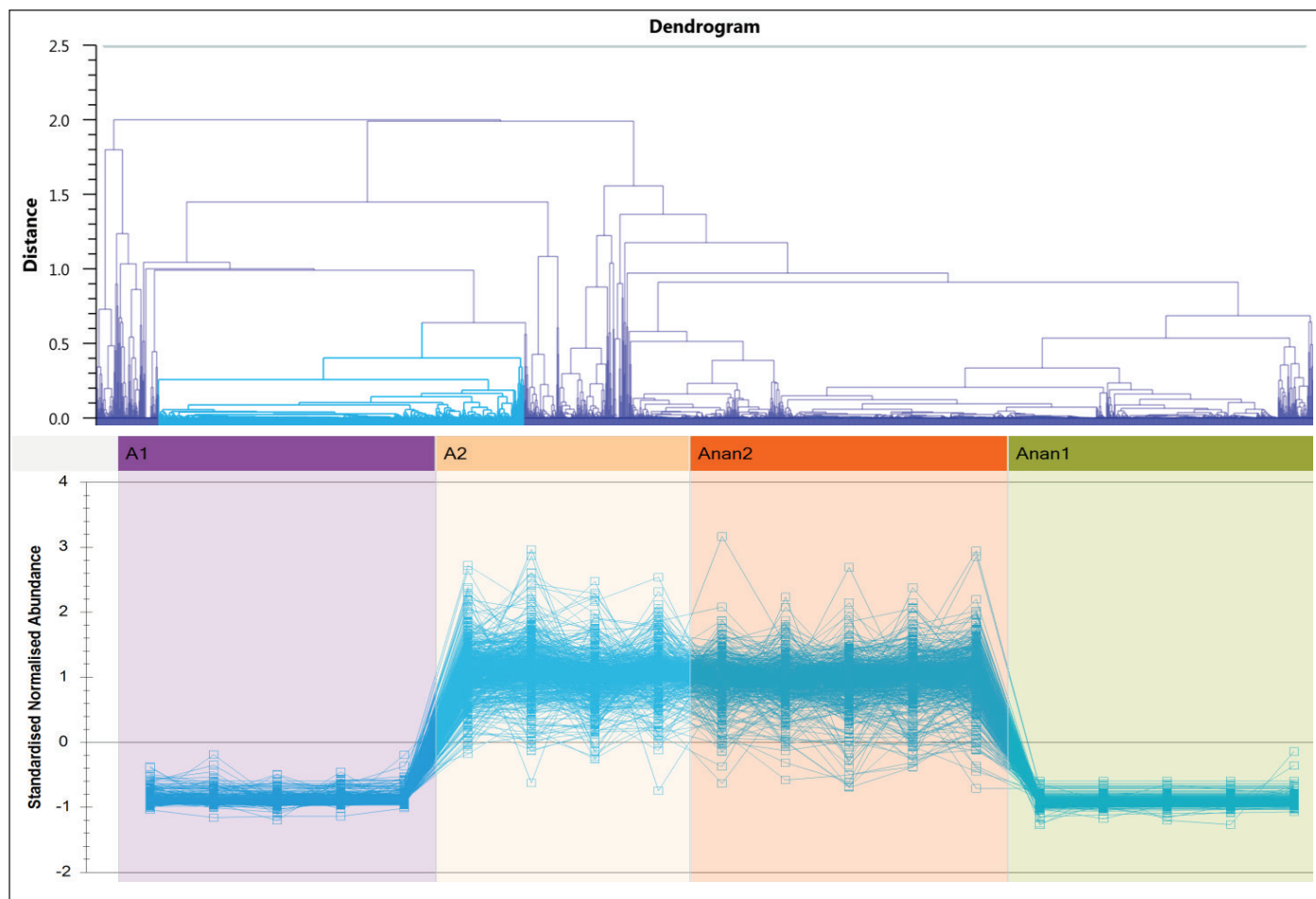


Figure 9. Expression and abundance profile of selected features according to their relative similarity between the different groups.

## CONCLUSIONS

The UPC<sup>2</sup>/MS FFA analysis described provides a simple and fast method with a significant reduction in analysis time compared to alternative techniques such as GC/MS, which requires FAME derivatization. In addition, the organic layer extract containing the lipids can be injected directly into the system, omitting the need for solvent exchange for compatibility with reversed-phase LC methods.

Saturated and unsaturated FFA containing C<sub>8</sub> to C<sub>36</sub> carbons were separated and determined, including low volatile very long chain fatty acids (>24 carbon atoms) that have challenged GC/MS even after FAME derivatization. Data analysis and FFA identification was facilitated using TransOmics for Metabolomics and Lipidomics Software that adopts an intuitive workflow approach to performing comparative ACQUITY UPC<sup>2</sup>/Xevo G2 QToF MS metabolomics and lipidomics data analysis.

## References

1. Fahy E *et al.* Update of the LIPID MAPS comprehensive classification system for lipids. *J Lipid Res.* 2009; 50: S9-S14.
2. Rustan AC, Drevon CA. *Fatty Acids: Structures and Properties.* eLS. 2005.
3. Allard B, Templier J. Comparison of neutral lipid profile of various trilaminar outer cell wall (TLS)-containing microalgae with emphasis on algaenan occurrence. *Phytochemistry.* 2000; 54: 369-380.
4. Stump CL, Goshawk J. The MarkerLynx Application Manager: Informatics for Mass Spectrometric Metabonomic Discovery. Waters Application Note 720001056en. 2004.



# Waters

**THE SCIENCE OF WHAT'S POSSIBLE.®**

Waters, ACQUITY UPC<sup>2</sup>, Xevo, UPC<sup>2</sup>, and The Science of What's Possible are registered trademarks of Waters Corporation. UltraPerformance Convergence Chromatography, and TransOmics are trademarks of Waters Corporation. All other trademarks are the property of their respective owners.

©2013 Waters Corporation. Produced in the U.S.A.  
July 2013 720004763EN AG-PDF

**Waters Corporation**  
34 Maple Street  
Milford, MA 01757 U.S.A.  
T: 1 508 478 2000  
F: 1 508 872 1990  
[www.waters.com](http://www.waters.com)